wau

KIO
BODY DEVICE
USER GUIDE

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RoHS

Table of Contents

Technical Specifications.....	20
Purpose.....	21
Device Components	21
Operating Procedure	22
Usage Method.....	23
Cleaning & Maintenance	27
Storage.....	27
Charging	27
Precautions	28
Warranty.....	32

Technical Specifications


Parameter	Value
Model	TBW-001
Material	ABS plastic
Dimensions (H×W×D)	111×91×28 mm (body) 111×91×44 mm (body with EMS attachment) 111×91×37 mm (body with RF&LED attachment)
Weight	128 g (body) 96 g (EMS attachment) 93 g (RF&LED attachment)
Power Source	Lithium-ion battery
Battery Capacity	1050 mAh
Operating Voltage	3,7 V
Power	10,6 W
Charging Time	3 hours
Operating Time (fully charged)	60 minutes
Smart Timer	5 minutes
Water Resistance Level	IPX6
Modes	EMS Mode RF&LED Mode

Purpose

The KIO body device with EMS, RF & LED functions by WAU is designed for body contouring in the areas of the abdomen, waist, hips, and buttocks. The massager addresses issues such as loss of skin firmness and tone, and helps combat sagging and cellulite using three of the most effective technologies.

EMS uses low-frequency electrical currents that cause muscle contractions, improve blood and lymph circulation, tighten the skin, and enhance muscle tone. RF waves release thermal energy, tighten the skin, and promote collagen regeneration, making the skin firmer and more elastic. Red LED light therapy targets fibroblasts, promoting skin rejuvenation and combating signs of aging.

Device Components

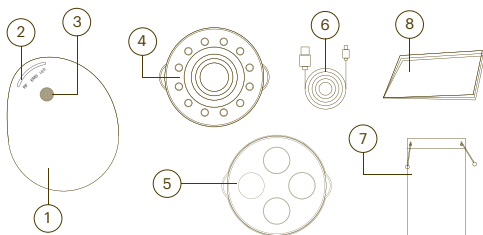1. Device body
2. Light indicator for intensity level/ battery level
3. Power button/ intensity adjustment
4. RF&LED attachment
5. EMS attachment

6. USB-C charging cable
7. Storage pouch
8. User guide


Package Contents

Device body x1, EMS attachment x1, RF&LED attachment x1, USB-Type C charging cable x1, storage pouch x1, user guide x1.

Operating Instructions

1. The KIO Body Device comes with two interchangeable heads: a large ball head for EMS mode and a ring head with small balls for RF & LED mode.

To install a head, align the groove on the back of the massage head with the raised section inside the attachment socket on the device body. Insert and secure the head in place.

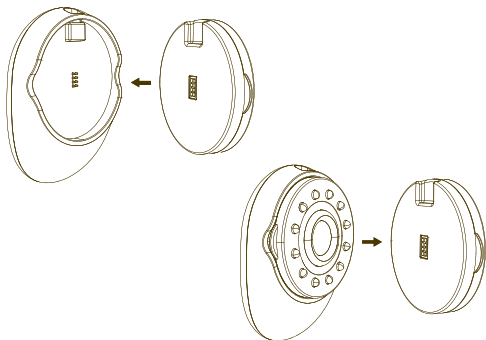2. Once the attachment is installed, turn on the device by pressing and holding the power button for 2 seconds. Depending on the installed head, the device will automatically activate either EMS or RF & LED mode. The intensity level indicator will light up.
3. Both EMS and RF & LED modes offer three intensity levels. To change the intensity, press the intensity adjustment button. The indicator will gradually fill as the intensity increases. After reaching the highest level, pressing the button again will cycle back to the lowest intensity. Note: In RF & LED mode, only the RF intensity changes. The LED intensity remains constant.
4. The device is equipped with a smart timer and will automatically power off after 5 minutes of operation.
5. To turn off the device, press and hold the power button for 2 seconds.
6. To remove the attachment, hold the device with one hand and firmly grip the tabs on the massage head with the other. Pull the attachment straight out with force to detach it from the device body.

Usage Method

General Recommendations

- Use the device only on clean skin. Make sure to remove any excess oils or impurities from the body before use.
- If your conductive gel or skincare product has been absorbed into the skin during the massage, reapply it to ensure smooth gliding of the massager and to enhance the treatment's effectiveness. As you move between different areas of the body, apply conductive gel to the specific area before starting the massage.

- The KIO device is equipped with a motion sensor. It is essential to keep the device moving during use. If left stationary on the skin for an extended period, the device will automatically stop. It will resume operation once movement is detected.
- Avoid repeated movements over the same area for an extended time, and do not keep the device stationary in one spot. Pain thresholds and skin sensitivity vary from person to person. Holding the device in one area without movement may cause minor burns. Do not stay on the same spot for longer than 2–3 seconds.
- Skin sensitivity and treatment effectiveness vary by individual and treatment area. Results and sensations may differ from person to person.

EMS Mode

1. Apply a conductive gel or your preferred water-based gel-like skincare product to the skin. This ensures optimal conductivity of the electrical currents and maximizes treatment results.
2. Attach the massage head suitable for EMS mode.
3. Turn on the device and, for first-time use, select the lowest EMS intensity level.

4. Waist Massage: Move the device along the sides of your waist, from the outer area toward the center, for approximately 5 minutes.

5. Leg Massage: Glide the device along the thighs or calves, moving upward in the direction shown. Continue for 5 minutes.

6. Arm Massage: Move the device from the elbow up toward the shoulder joint for 5 minutes.

7. Abdomen Massage: Move the device in a zigzag motion vertically up and down, avoiding the navel area. Continue for 5 minutes.

8. While massaging, apply firm pressure to ensure full contact between the electrode surface and the skin.

9. After completing the session, turn off the device, remove and clean the massage head. Leave the conductive gel or skincare product on the skin until fully absorbed. Then cleanse the skin if desired.

RF Mode

1. Apply a conductive gel or your water-based skincare product with a gel-like texture to the skin. This will ensure better current conductivity for optimal results.

2. Attach the appropriate head according to the mode you want to use.

3. Turn on the device. For the first use, select the initial RF&LED intensity level.

4. Choose the target area and perform the massage movements listed below.

5. Using gliding motions, move along the texture of the skin for 5 minutes.

6. Glide horizontally for 5 minutes.

7. Move in circular motions clockwise for 5 minutes.

8. Repeat circular motions counterclockwise for 5 minutes.

[Note]: The core function of RF is to heat the skin, stimulating collagen regeneration and subcutaneous fat metabolism. If the heating feels insufficient, reduce the gliding speed. If the temperature feels too high, increase the speed to maintain a comfortable sensation.

9. After the procedure, turn off the device, remove and clean the attachment. Leave the conductive gel or skincare product on the skin until fully absorbed, then cleanse the skin.

Frequency and Duration of Use

- Use one or both modes depending on your needs.
- The device is equipped with a smart timer: it will automatically turn off after 5 minutes. This marks a single session segment. If you wish to continue the massage, simply restart the device.
- When treating the same area, allow at least 1 day of rest between sessions.
- Use the device 3 times a week for at least 6–8 weeks. After this period, we recommend continued regular use to achieve the best results. Frequency and duration of use can be adjusted according to your individual needs.

Cleaning & Maintenance

- Make sure the device is turned off before cleaning.
- After each use, clean the attachments and the body of the massager. This device has an IPX6 waterproof rating, which allows rinsing it under running water. However, do not submerge the device in water or any other liquid.
- After cleaning, wipe off any remaining water with a dry cloth.
- Do not use metal scrubbers, abrasive cleaners, chemicals, or corrosive liquids to clean the device.

Storage

- Store the device in a dry place, protected from dust.
- Keep the device away from high temperatures.
- Store the device out of reach of children.
- If you plan not to use the device for an extended period, charge it before storing.

Charging

- When the battery level drops below 10%, the battery indicator starts flashing rapidly.
- At 0% battery, the indicator flashes 5 times before the device powers off.
- Connect the charging cable with a Type-C connector to the device's port. Plug the USB end of the cable into a power adapter.
- During charging, the battery indicator gradually fills and then turns off.
- The device does not operate while charging.
- When the battery is fully charged (100%), the indicator remains fully lit.

Safety Precautions

Please read this user guide carefully before use and retain it for future reference.

Warnings

- This device is not a toy. Keep it out of reach of children.
- Use the device with care and only as instructed. Failure to follow the guidelines may result in unwanted effects.
- Prolonged, unsupervised use of the device on the same area of the body may cause potential skin damage.
- Do not use the device in combination with flammable cosmetic products.
- Do not use the device around the eye area.
- Inspect the device for any damage before each use. Do not use the device if it is malfunctioning, broken, or damaged.
- Perform a sensitivity test 24 hours prior to first use. Apply a small amount of contact gel to the inner forearm to ensure neither the skincare product nor the device causes a skin reaction.
- Do not use the device for extended periods on body areas where you do not wish to reduce volume.
- Using the device in dry environments may cause static electricity. This is a normal occurrence and does not affect the functionality of the massager.
- Keep the USB-C charging cable and the charging port dry. Do not operate the device while charging.
- Use caution when applying contact gel or skincare products. Oils and acids in cosmetic formulas may damage the device's surface.
- Consult your physician before use.
- Discontinue use if you experience any discomfort, allergic reactions, skin irritation, dermatitis, or any other side effects.

Contraindications

Do not use this device if you:

- Have a pacemaker or any other implanted electronic medical device;
- Suffer from allergies;
- Have diabetes, lupus, porphyria, congestive heart failure, blood clotting disorders, or autoimmune diseases (e. g., HIV or AIDS);
- Have active implants such as a pacemaker or internal defibrillator;
- Have arrhythmia or any other heart conditions;
- Have permanent metal implants or body piercings in the area of device use;
- Have any unidentified skin lesions;
- Are undergoing radiation therapy or chemotherapy;
- Have conditions caused or aggravated by heat exposure.
- This device is suitable for use on stretch marks caused by obesity or pregnancy; however, it is not intended for scar treatment.

Do Not Use This Device If You Are Taking Any of the Following Medications:

- If your skin has been treated recently or currently with alpha hydroxy acids (AHA), beta hydroxy acids (BHA), topical isotretinoin, or azelaic acid;
- If you have taken any form of Roaccutane (isotretinoin) within the past 6 months;
- If you are taking photosensitizing medications or any other prescribed drugs. Please carefully review the instructions for your medication before using the massager. If the medication's instructions indicate the possibility of photoallergic or phototoxic reactions, or advise avoiding sun exposure, do not use this device;
- If you are taking anticoagulant medications, including high doses of aspirin, use the device only one week after completing treatment.

Additional Information:

— This device is intended for use on the hands, waist, abdomen, thighs, and buttocks only. Do not use the device on the face, neck, chest, genital area, or on any body areas that have undergone plastic surgery within the last 3 months.

— Do not use the device in LED therapy mode if your skin is sensitive to light or prone to rashes or allergic reactions.

— Do not use the device in LED therapy mode if there are infections, eczema, burns, folliculitis, open wounds, bruises, cold sores, skin lesions, or hematomas in the area of application.

— Do not direct the light source into the eyes or look directly at the light emitted by this device.

— Do not attempt to disassemble, repair, or modify the device. Use this device only for the purposes specified in this user guide.

Malfunctions and ways to eliminate them

Malfunction	Reasons
The device does not work after being turned on.	Low battery level or completely discharged battery.
How to fix	
Charge the device according to the instructions.	

Malfunction	Reasons
The heating is not felt	1. You haven't applied enough contact gel. 2. You're moving the device across your body too quickly. 3. The device has stopped working.
How to fix	
1. Apply the conductor gel or your skin care product to the skin. 2. Move the device slowly over the body. 3. Continue massaging at a moderate speed.	

Recommendation: If the device does not work properly despite the measures you have taken, please stop using the massager and consult with an organization authorized by the manufacturer to accept consumer claims.

Warranty

The KIO Body Device by WAU is covered by a 1-year warranty. During the warranty period, if technical defects are found, the product is subject to replacement or repair.

This warranty does not cover malfunctions caused by improper use, transportation damage, or products that have been repaired by unauthorized persons.

KIO

BODY DEVICE

Model: **TBW-001**
Material: **ABS Plastic**
Dimensions (HxWxD): **107 x 88 x 32 mm**
Weight: **317 g**
Voltage: **3.7 V**
Power: **10.6 W**
Made in: **PRC**

Authorized manufacturer representative:
I. E. Umantsev A. N. Leningradsky Prospekt 37, Office
14-14, Moscow, 125167, Russia
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS